aquapick[®]

PORTABLE ORAL IRRIGATOR

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Advanced Oral Irrigation System
Sistem avansat de irigare orală

IMPORTANT SAFEGUARDS

Please read the safety instructions before using this product. The safety instructions provided here help you to use this product correctly, and also to prevent you from injuring yourself or others. When using electrical products, especially when children are present, the basic safety precautions should always be followed.

READ ALL INSTRUCTIONS BEFORE USE

1. DANGERS

- ❶ Toddlers, infants, Children should not be allowed to use this device. Please keep it out of reach of children.
- ❷ Do not use if the power cord or plug was damaged.
- ❸ Always unplug the device after using.
- ❹ Do not touch plug with wet hand.
- ❺ Do not place the device in water or other liquids.
- ❻ Do not place the device where it may fall or be pulled into a tub or sink.
- ❼ Do not use the device while bathing.

2. WARNINGS

- ❶ The device can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under supervision.
- ❷ Children should not play with the appliance.
- ❸ Use only water or neutral cleaning agent to clean the device. Never use caustic or abrasive agents (e.g., vinegar or lime scale removers), since they may damage the device.
- ❹ Keep jet tip away from small children. The jet tip is small enough to be swallowed and can become lodged in the throat.
- ❺ Do not direct water under the tongue, into the ear, nose or other delicate areas.
- ❻ Do not use for any other purpose except oral cleaning.
- ❼ Do not fill the tank with warm water over than 40 °C.
- ❽ People who experience tooth or gum pain may not be able to use this device and should consult their dentist.
- ❾ Do not press the jet tip too close against your teeth or gums, doing so may injure your gums.
- ❿ Please fill the tank with pure water.

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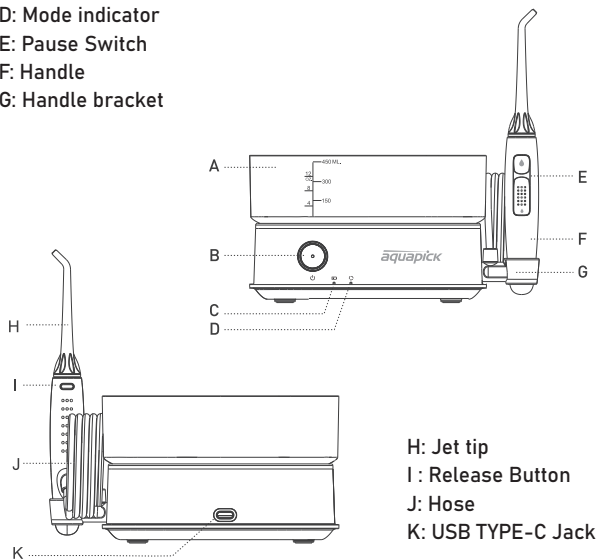
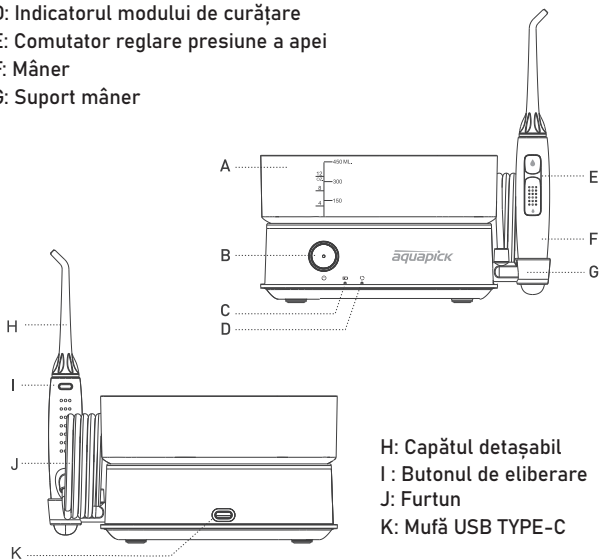
HOW TO USE

Getting Started

Place the device on a flat surface without any sharp objects and make sure the cord is dry and has no damage. The device is powered by internal battery and charged with USB TYPE-C jack connector.

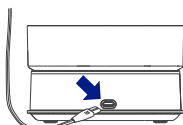
Fig.1

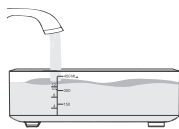
Fig.2-1

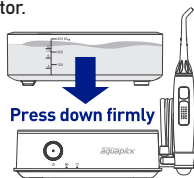
Fig.2-2

Charging the device

- 1 The machine must be charged for first time.
- 2 Insert the plug into the USB TYPE-C jack for charging to make sure the Charge indicator flashing until the charge indicator turn to white then unplug. (Fig.1)
 - ※ A full charge requires 4 hours, don't charge too long.
 - ※ If the machine stops working due to excessive heat during use, please wait 2 hours before use.

Using the device

- 1 Take the Water Tank off to fill with proper liquid, no warmer than 40 °C, then replace the tank onto the top of main unit and **press down firmly**. (Fig. 2-1 & Fig. 2-2)
 - ※ **Please note that if the main body and the water tank do not join properly, it may malfunction.**

Inserting& Removing Jet tip

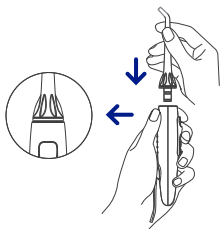
Fig.3

Inserting jet tip

Insert tip into the center of the handle, audibly indicating when the jet tip attached to the handle is in place. Also try to pull the tip out gently to confirm the tip is locked firmly. (Fig.3)

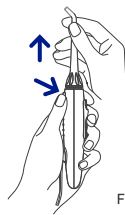
Fig.4

Removing jet tip

Press the Release Button on the top end of the handle, then pull the jet tip out lightly by hand. (Fig.4)

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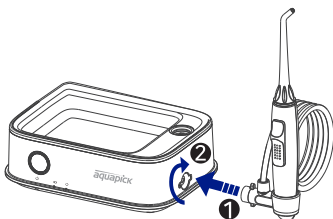
Fig.5

Inserting Handle

Insert handle bracket into the bottom of the machine, and clockwise rotation gently to confirm the handle is locked firmly. (Fig.5)

Removing Handle

Contrarotate the Handle bracket to unlock.

HOW TO USE THE ORAL IRRIGATOR

Notes:

Do not power on the device before placing the jet tip in your mouth. First use, we suggest you choose the lower pressure, gradually increase the water pressure over time to the setting you need.

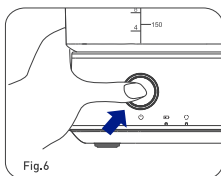
Fig.6

- 1 Hold the handle and point the jet tip to sink, then press down the Power Button to power the device on. (Fig.6)

Press the Power button again within 10 seconds to select the mode you desired. (Clean or Massage), default Clean mode for the first time.

Working mode: Clean mode/ Massage mode

Mode indicator:

Blue light is on when working at Clean mode.

Green light is on when working at Massage mode.

Power/Charge indicator:

White light is flashing during charging.

Red light is flashing when battery is low.

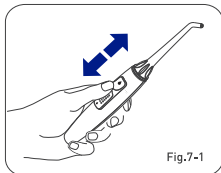
Fig.7-1

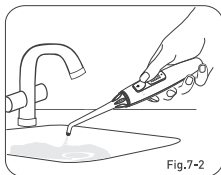
Fig.7-2

- 2 Push the Pause Switch to adjust water pressure, divided into 5 pressure gears. (Fig.7-1)
- 3 Adjust the direction of jet tip by rotating the jet tip kob. (Fig.7-2)

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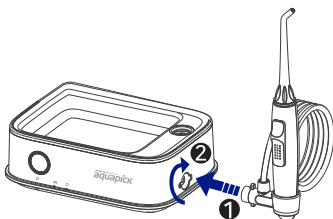
Fig.5

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
- ④ Bend over to the sink and place the tip in your mouth, and then push the Pause Switch upwards slowly to enjoy irrigating. (Slightly close your mouth to avoid the water splashing.) (Fig.8)


- ⑤ Position the water jet flow at 90-degree angle to the teeth and gums. (Fig. 9)
- ⑥ Move the water jet slowly along the teeth. The unit will continue to work until 2-minute cycle is complete or you turn it off.

AFTER USE

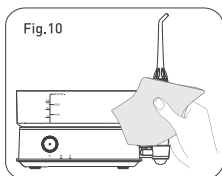- ① Take off the water tank and dispose of water in it.
- ② Push the Pause Switch upwards, and power on to drain the water inside the device, then push the Pause Switch downwards.

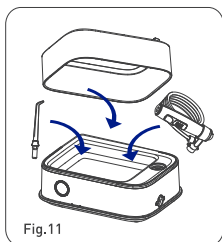- ③ Turn off the device, rinse the water tank and wipe it with a soft, dry cloth.
- ④ Press down the Release Button to remove the jet tip and then place the jet tip and the handle into the main body and cover it with a water tank. (Fig.11)

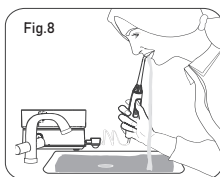
Fig.8

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TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
Unit does not work	Motor not running.	Make sure that the charging is done well.
Water pressure is very low	Dirt stuck in the jet tip.	Clean or replace with a new jet tip.
	The jet tip is deformed.	Replace with a new jet tip.
	Dirt stuck in the tank valve.	Clean the valve and fill the tank with pure water.
Water does not come out	The water bottle is not well bonded to the main body.	Insert the water tank into the top of the body and press firmly to ensure tight engagement.
	Water pressure under the lowest setting.	Slide the switch to adjust pressure.
	Tank valve not open.	Press down the tank, or reassemble the tank valve.
Tank leaks	The valve was lost.	Replace with a new tank valve.
	The valve isn't install in place.	Take off the valve and reinstall again.

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